


# LeoMed

## Ligation Device


UNIQUE  
DESIGN

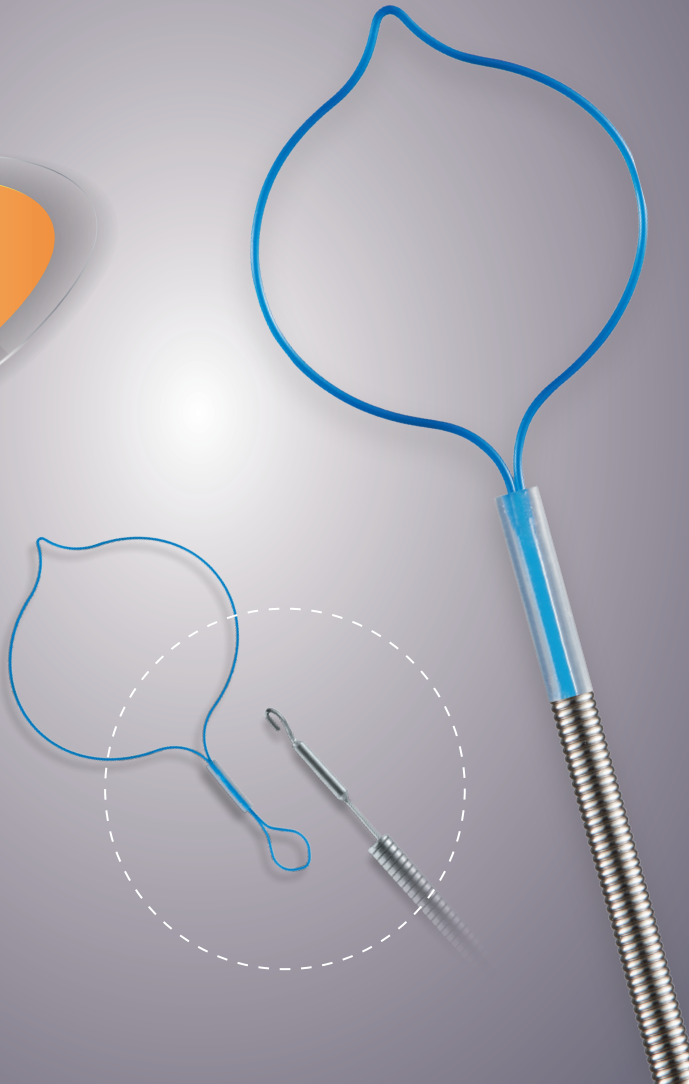

### Intended Use:

Used for endoscopic ligation of polyps.

### Product Feature:

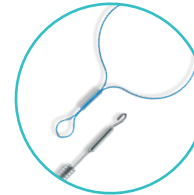

#### Split type of delivery system

With split design, the delivery system and the loop do not need to be pre-assembled, and multiple ligation loops can be operated inside human body.

#### Various Loop sizes

Available in various clinical sizes.



#### Unique double circle design

The double circle structure carefully designed according to clinical needs is more conducive to purse-string suture with a single-channel endoscope, which is simple and convenient to use.

#### Excellent material selection

Imported nylon loop provides moderate hardness to ensure safe and effective ligation operation.

### Specification Suggestion

If you need other models, please contact our sales team.

#### ▲ Split Type

Loop

REF	Open Diam.(mm)	Height H (mm)	Silicon Tubing Diam. D (mm)
Loop15B0	15	36	2.1
Loop20B0	20	40	2.1
Loop30B0	30	55	2.1
Loop40B0	40	70	2.1

#### ▲ Preassembled(Disposable)

REF	Open Diam.(mm)	Effective Length (mm)
Loop15B0-LD230D0	15	2300
Loop20B0-LD230D0	20	2300
Loop30B0-LD230D0	30	2300
Loop40B0-LD230D0	40	2300